

LUTRONIC



INFINI™
HIGH INTENSITY FOCUSED RF
The New Fractional


Precisely Delivered Radiofrequency
Through Insulated Microneedles for
Consistent Outcomes You and Your
Patients Will Love



INFINI RADIOFREQUENCY MICRONEEDLING

INFINI is unique technology that puts you in control of your procedure. You control the treatment depth and energy levels to enable a multi-layered approach for customized treatments, regardless of skin type or tan. Now you can provide the best solution that will keep your patients returning and referring you to new patients.

EACH TREATMENT SESSION PLACES 50,000+ MICRO-COAGULATION ZONES.

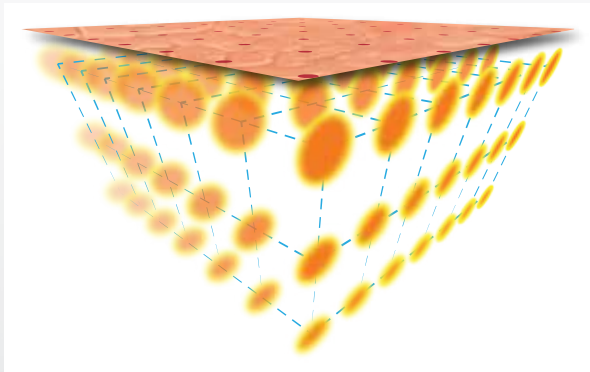

WHY INFINI WITH RF ENERGY?

- Creates thousands of coagulation zones
- Precise energy delivery
- Automatic settings
- Safe on all skin types
- Low cost consumables
- Minimal downtime (makeup after 24 hours)
- Easy to delegate
- 34 gauge microneedles for patient comfort


PROVEN RESULTS FOR ACNE SCARS, WRINKLES, TEXTURE, LAXITY¹



Courtesy of Steve Weiner, MD




Courtesy of Edmond Griffin, MD

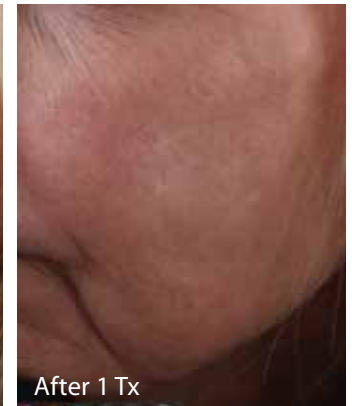

Courtesy of Omar Ibrahim, MD, PhD




Courtesy of Jason Pozner, MD

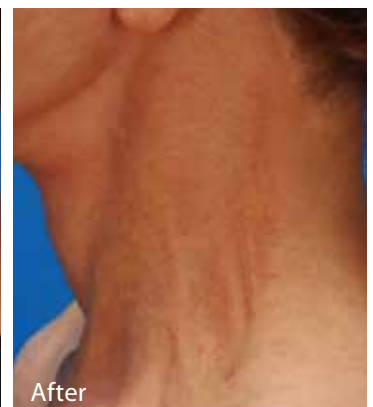


TECHNOLOGY BREAKTHROUGH OF THE DECADE

MICRONEEDLE RF TREATMENT PRODUCES WELL DEFINED ZONES OF COAGULATIVE DAMAGE

Infini is clinically proven, with studies to back the technology. The histology images below show the discrete coagulative injury zones in the dermis while sparing the epidermis. The natural healing process replaces the damaged zones with new collagen and elastin.

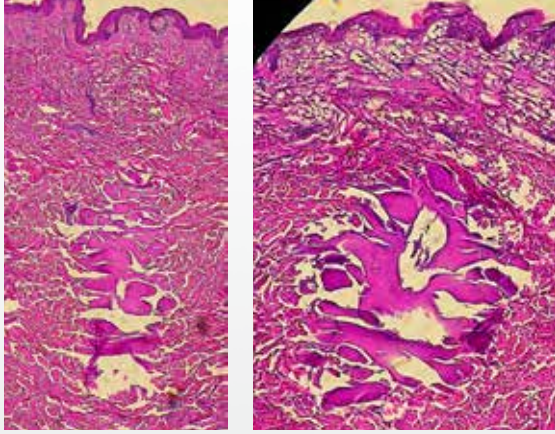

Ex vivo human tissue shows epidermis is left intact, while variable depth coagulation zones are present.

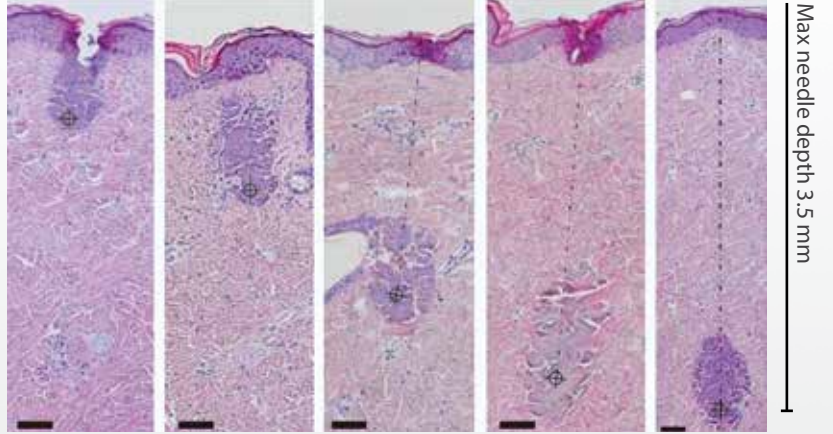

0.5 mm 1.0 mm 1.5 mm 2.0 mm 3.5 mm


All depths treated at Level 4 for 100 ms

Courtesy of SB Cho, MD, Dermatologist, Yonsei University, S. Korea


IDEAL FOR THE NECK¹



Courtesy of Oscar Hevia, MD

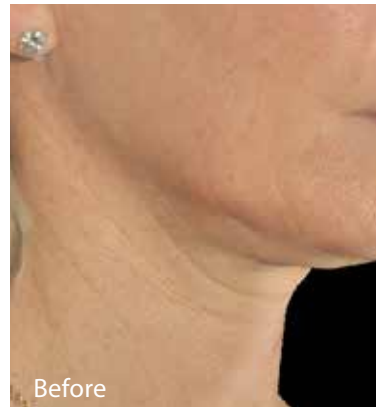


Courtesy of Steve Weiner, MD

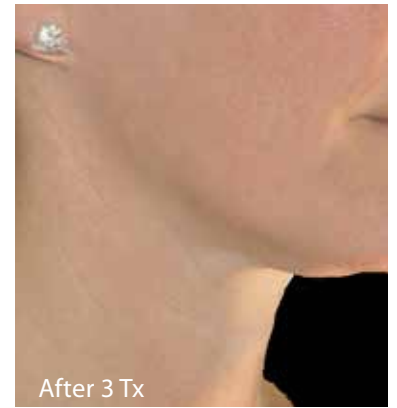

Courtesy of Sarah Mess, MD



Courtesy of Matteo T. Clementoni MD

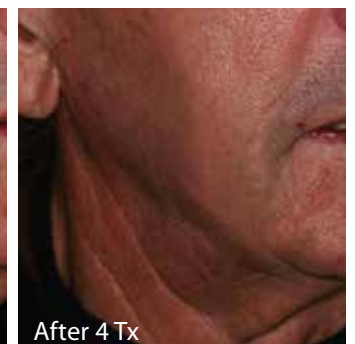

LUTRONIC

¹Infini is FDA cleared for the treatment of facial wrinkles

DELIVERING OPTIMIZED OUTCOMES IS AS EASY AS 1, 2, 3¹

1

SELECT TREATMENT ZONE


2

SELECT PASS DEPTH

3

SELECT INTENSITY

¹Manual Mode still available

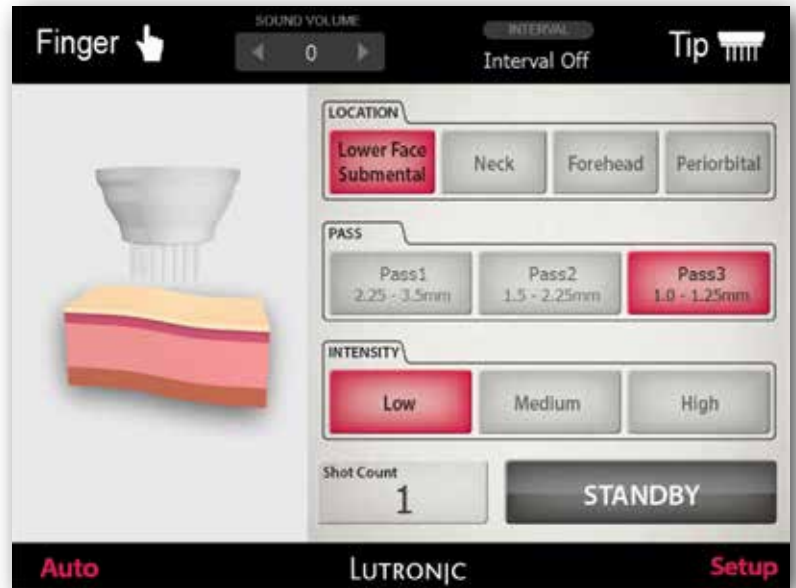

Proprietary optimized pin to pin RF transmission across 7X7 array of 49 insulated microneedles



WHAT PHYSICIANS ARE SAYING



"I've been involved in the laser device business for many years. When it really came down to the science, I believe INFINI provides the most tightening for the effort and least amount of downtime."

Jason Pozner, MD



"I've been using INFINI for about a year, we have been really excited about the results we have seen for acne scars and neck* rejuvenation. I think one of the things that stands out with INFINI as opposed to other devices is minimal downtime. You can finish a treatment and just be a little pink the next day."*

Joel Cohen, MD

*Infini is FDA cleared for the treatment of facial wrinkles

SYSTEM SPECIFICATIONS*

System Type	Bipolar RF (Radiofrequency)
Weight	28 kg
Dimension	362 mm (W) x 409 mm (L) x 1713 mm (H)

lutronic.com | 888-588-7644

©2016 LUTRONIC Group of Companies. All rights reserved. LUTRONIC, its logo, and INFINI are trademarks or registered trademarks of the LUTRONIC Group of Companies. *Specifications are subject to change without notice. This product or the use of this product is covered by one or more US and foreign patents or patent application pending. INFINI_4100149300_US_v5

LUTRONIC